

ProClense plus SCD™

Super Concentrated Alkaline Instrument Detergent

Description:

ProClense plus SCD is a super concentrated alkaline detergent used for cleaning instruments and devices in automated equipment. It is super concentrated to reduce package weight and storage space.

Super Concentrated:

- Clinically tested for effective cleaning at ultra low dose settings.
- Reduces weight of containers and storage space.
- Delivers critical cleaning performance with reduced cost per cycle.

Alkaline Detergent:

- Indicated for use when an alkaline detergent is preferred.
- Powerful cleaning action for greasy soils and contaminated lab ware.
- Tested for compatibility with stainless steel, titanium, gold, brass, copper, heat stable plastics and rubber goods.


Enhanced Cleaning Action in Automated Equipment:

- Non-foaming for easy rinsing.
- Formula promotes ideal pump pressure and spin action.
- Chelating agents for best performance over wide range of water quality and operating conditions.



More Power
with Less Cost & Weight!

2.5 gallon drums are equivalent to 30 gallon drums of regular dosing detergent.

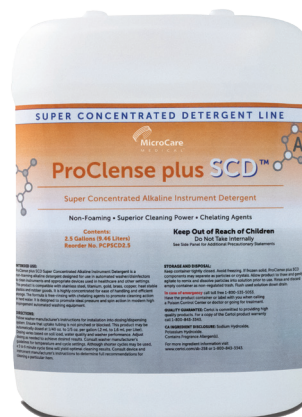

Chelating agents improve performance across a wide range of water quality.

Product Specifications

Dilution Ratio	1/40 to 1/5 oz. per gallon (0.2 mL to 1.6 mL per liter)
pH	13.5 to 14.5 (in concentrate)
Temperature of Use	150°F (66°C) or according to washer manufacturer's settings and specifications
Odor	Non-specific
Color	Natural clear yellow
Shelf Life	3 years from date of manufacture
Packaging	2.5 gal. drum, 1 each

Description	Case Qty	Order#
2.5 Gallon Drum	1	PCPSCD2.5-1

Product Capabilities:

ProClense plus SCD contains surfactants and other components to rapidly dissolve a wide range of surgical soils. The alkaline formula is highly effective on fats, mucous and blood. It is appropriate for use on heat tolerant instruments and devices indicated for cleaning in automated washing equipment.

Intended Use:

ProClense plus SCD Super Concentrated Alkaline Detergent is non-foaming and formulated to promote ideal pressure and spin action in automated washer/disinfectors. This product is used to clean instruments and appropriate devices used in healthcare and other settings. The alkaline formula is anti-corrosive and tested for material compatibility with stainless steel, titanium, brass, gold and copper. All applications should be according to manufacturer's instructions.

Precautions:

Use Certol SCD super concentrates in cleaning equipment designed to dose at more precise volumes from 1/5 oz. per gal. to 1/20 oz. per gal. (0.2% – 0.05%).

Availability:

Stocked by MicroCare Medical distributors. Contact MicroCare Medical for a distributor near you or visit MicroCare.com/Certol.

MicroCare Medical
6120 East 58th Avenue
Commerce City, CO 80022
(800) 843-3343



Visit Us:
MicroCare.com/Certol